


三七 (中国药典)

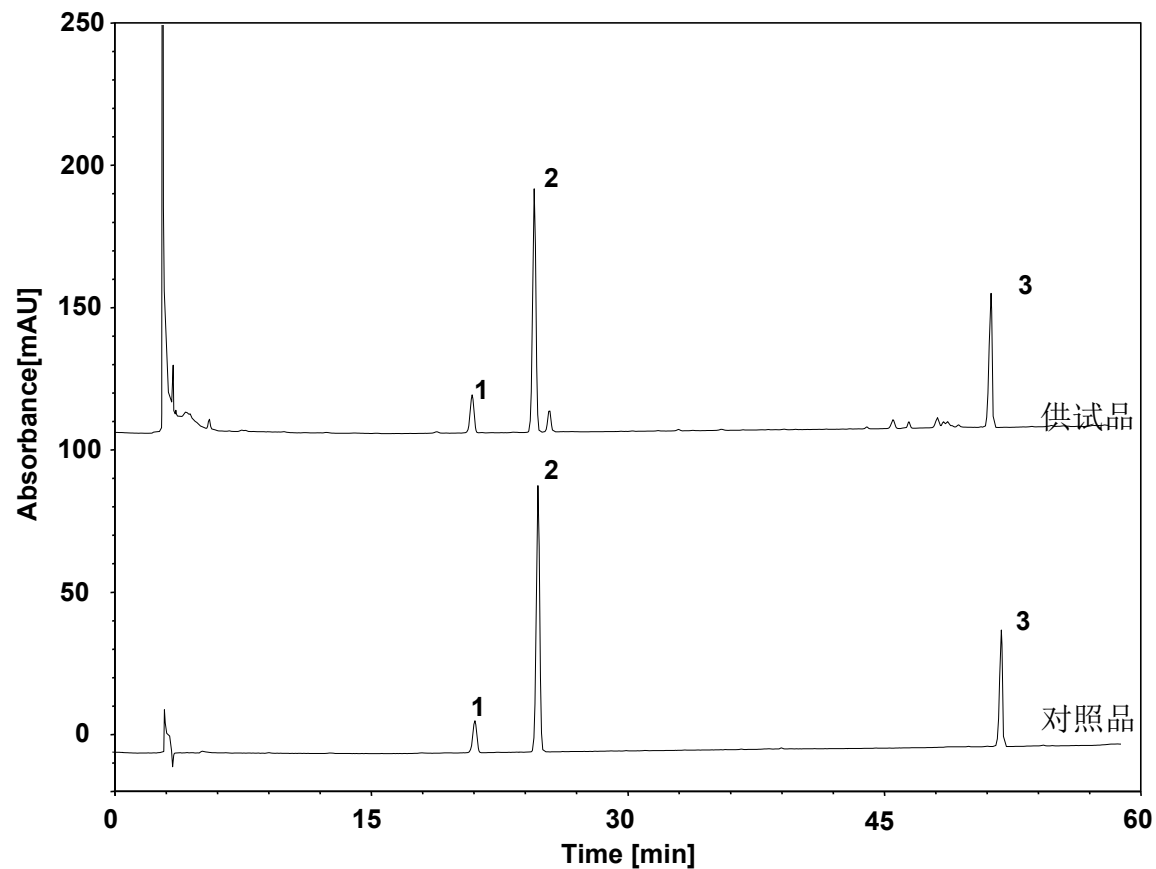

Column: ChromCore 300 C18, 5 μ m

Dimension: 4.6 \times 250 mm

Mobile Phase: A) 乙腈

B) 水

Gradient: t (min) %A %B

0 19 81

12 19 81

60 36 64

Flow Rate: 1.0 mL/min

Temperature: 30 $^{\circ}$ C

Injection: 10 μ L

Detection: UV 203 nm

Sample: 三七

Peaks:

1. 三七皂苷 (三七皂苷R1)

2. 人参皂苷 (人参皂苷Rg1)

3. 人参皂苷 (人参皂苷Rb1)