



# 3- (3-三氟甲基苯基)丙酸

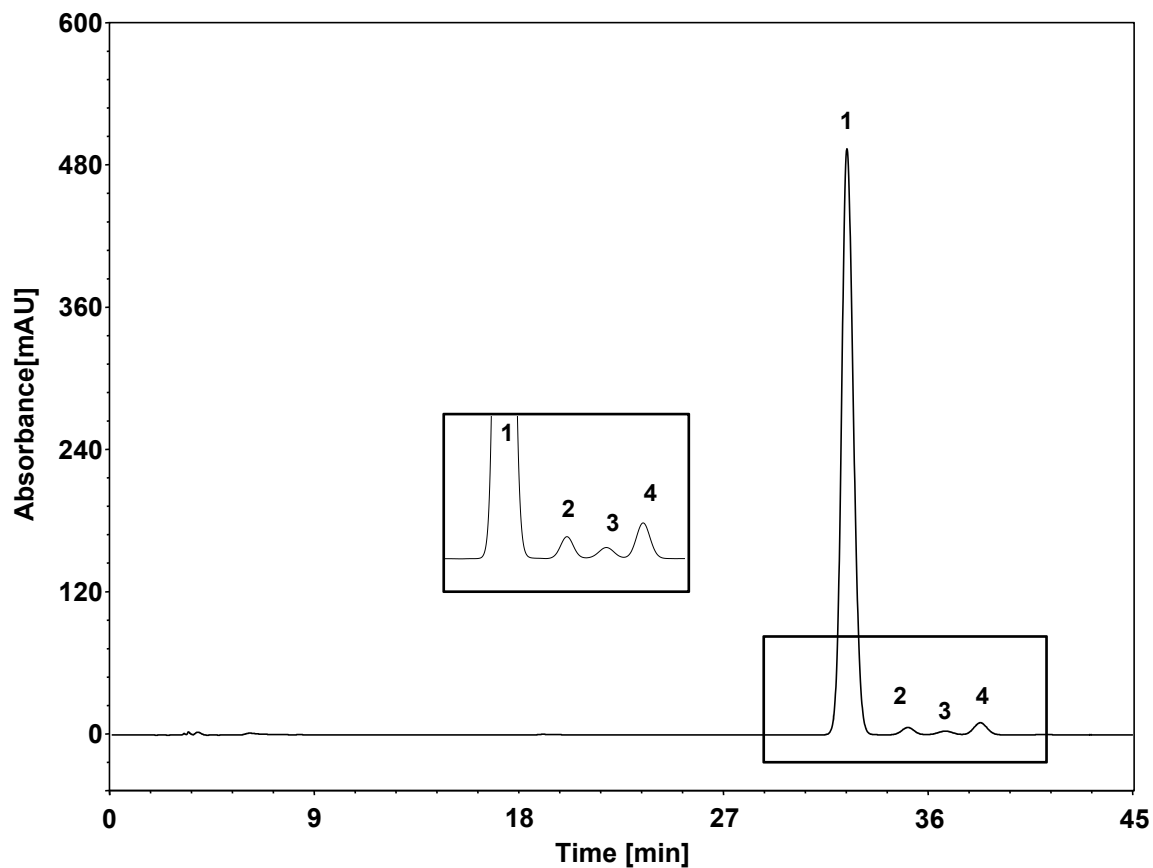


Column: ChromCore PFP, 5  $\mu$ m  
Dimension: 4.6  $\times$  250 mm  
Mobile Phase: 50/50 v/v 甲醇/0.1%磷酸水  
Flow Rate: 1.0 mL/min  
Temperature: 40  $^{\circ}$ C  
Injection: 10  $\mu$ L  
Detection: UV 210 nm  
Sample: 3- (3-三氟甲基苯基)丙酸混标  
Peaks:


1. 3- (3-三氟甲基苯基)丙酸
2. 杂质 D
3. 杂质 G
4. 杂质 A

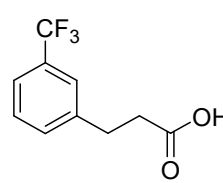


3- (3-三氟甲基苯基)丙酸

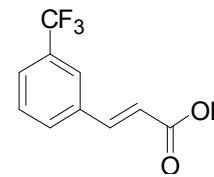

杂质A

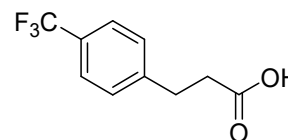

杂质D

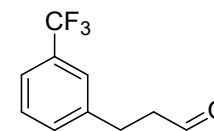

杂质G