


醋酸氟轻松乳膏 (中国药典)

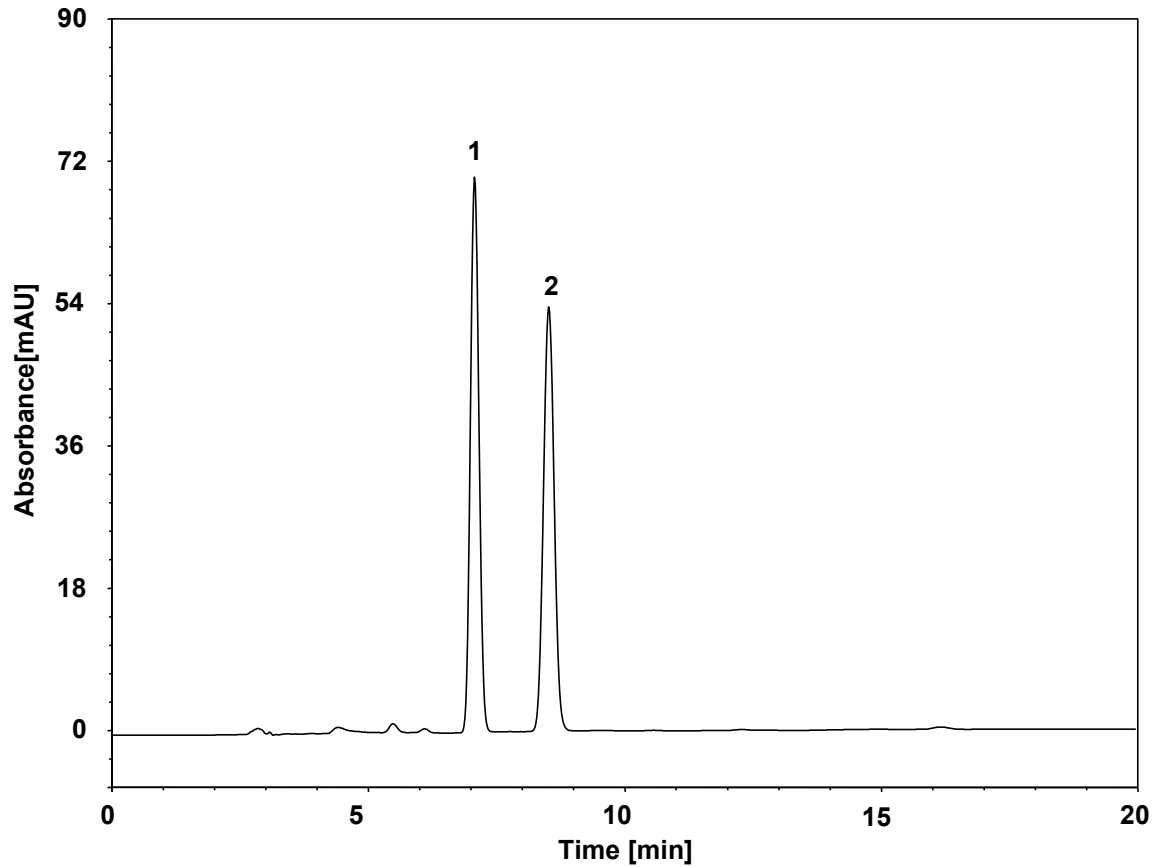

Column: ChromCore C18, 5 μ m
Dimension: 4.6 \times 250 mm
Mobile Phase: 60/10/30 v/v/v 甲醇/乙腈/水
Flow Rate: 1.0 mL/min
Temperature: 25 $^{\circ}$ C
Injection: 20 μ L
Detection: UV 240 nm
Sample: 醋酸氟轻松乳膏系统适用性溶液
Peaks:
1. 炔诺酮
2. 氟轻松