



炔丙基半胱氨酸

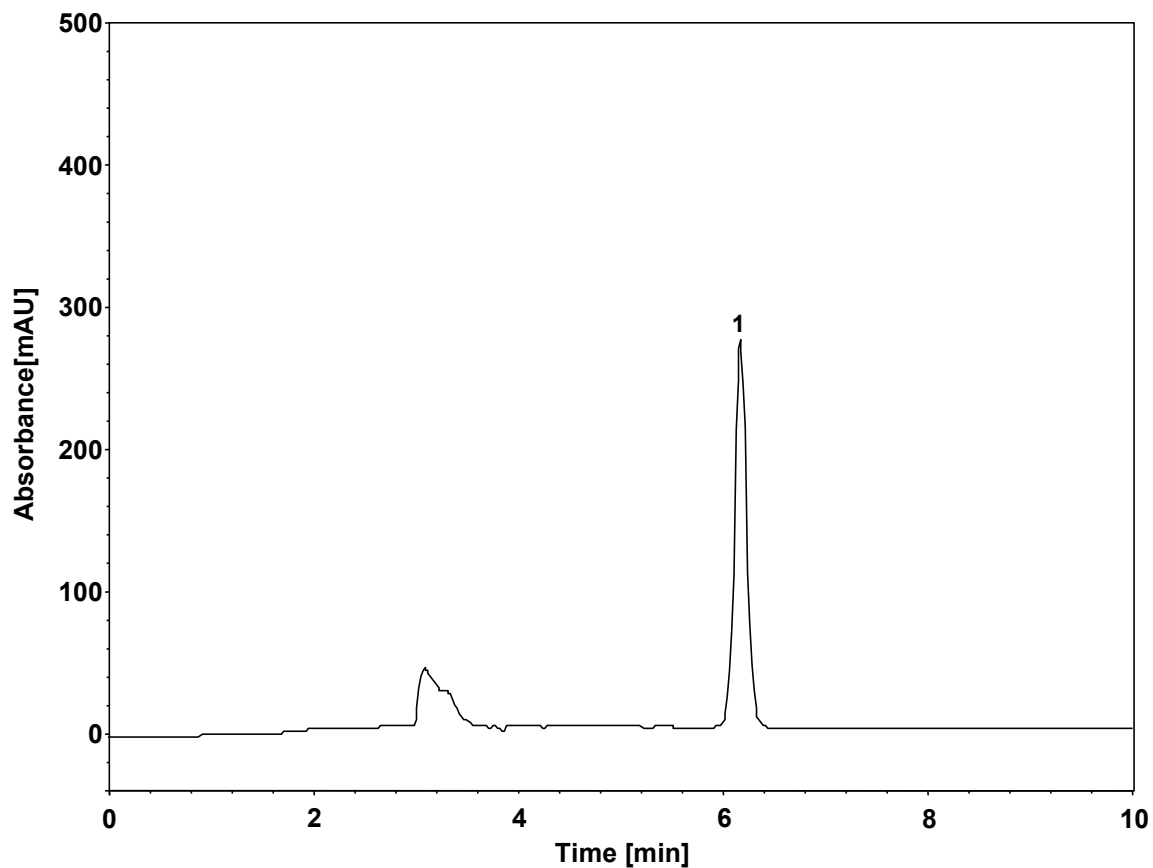

Column: ChromCore HILIC-Amide, 5 μ m
Dimension: 4.6 \times 250 mm
Mobile Phase: 70/30 v/v 乙腈/25 mM磷酸二氢钠溶液, pH6.0
Flow Rate: 1.0 mL/min
Temperature: 25 $^{\circ}$ C
Injection: 10 μ L
Detection: UV 210 nm
Sample: 炔丙基半胱氨酸 (0.5 mg/mL 50%乙腈水溶解)
Peak: 1. 炔丙基半胱氨酸

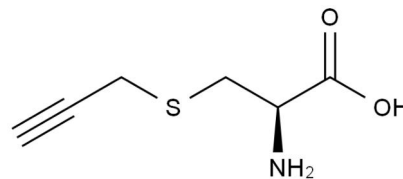

炔丙基半胱氨酸