




苯胺类在NanoChrom BP-5MS上的分离与检测 (EPA 8270)

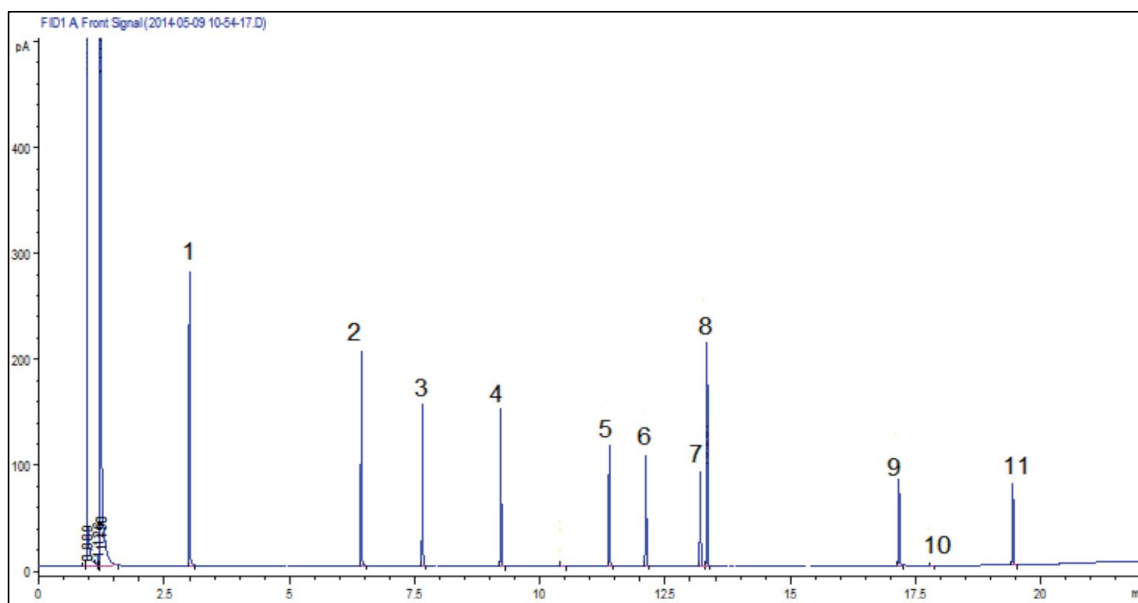

Column: NanoChrom BP-5MS
 Dimension: 30m×0.25mm×0.25μm
 Carrier: He, pressure: 12.2 psi
 Split: 60 mL/min
 Oven: 60 °C to 240 °C (hold 2 min) at 10 °C/min
 Inject: 0.1 μL
 Peaks:

1. N-亚硝基二甲胺
2. 吡啶
3. 苯胺
4. N-亚硝基二丙胺
5. 4-氯苯胺
6. 2-硝基苯胺
7. 3-硝基苯胺
8. 4-硝基苯胺
9. 二苯胺
10. 联苯胺
11. 3,3'-二氯联苯胺