




PBDE在NanoChrom BP-5HT上的分离与检测

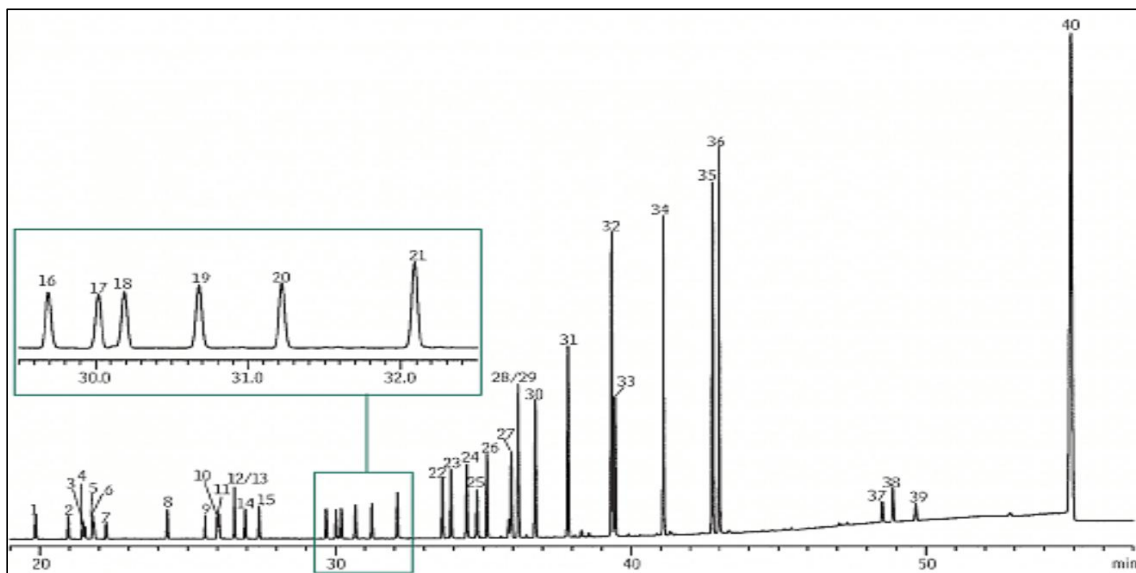

Column: NanoChrom BP-5HT
 Dimension: 15m×0.25mm×0.10μm
 Carrier: He, pressure: 12.2 psi
 Split: 60 mL/min
 Oven: 100 °C (hold 3 min) to 320 °C (hold 2 min) at 5 °C/min to 380 °C (hold 2 min) at 5 °C/min
 Inject: 1 μL

Peaks:	1. BDE-10	21. BDE-77
	2. BDE-7	22. BDE-100
	3. BDE-8	23. BDE-119
	4. BDE-11	24. BDE-99
	5. BDE-12	25. BDE-116
	6. BDE-13	26. BDE-118
	7. BDE-15	27. BDE-85
	8. BDE-30	28. BDE-155
	9. BDE-32	29. BDE-126
	10. BDE-17	30. BDE-154
	11. BDE-25	31. BDE-153
	12. BDE-28	32. BDE-138
	13. BDE-33	33. BDE-166
	14. BDE-35	34. BDE-183
	15. BDE-37	35. BDE-181
	16. BDE-75	36. BDE-190
	17. BDE-49	37. BDE-208
	18. BDE-71	38. BDE-207
	19. BDE-47	39. BDE-206
	20. BDE-66	40. BDE-209