


# 花粉祛痒止痛酊 (ZGB2019-31)

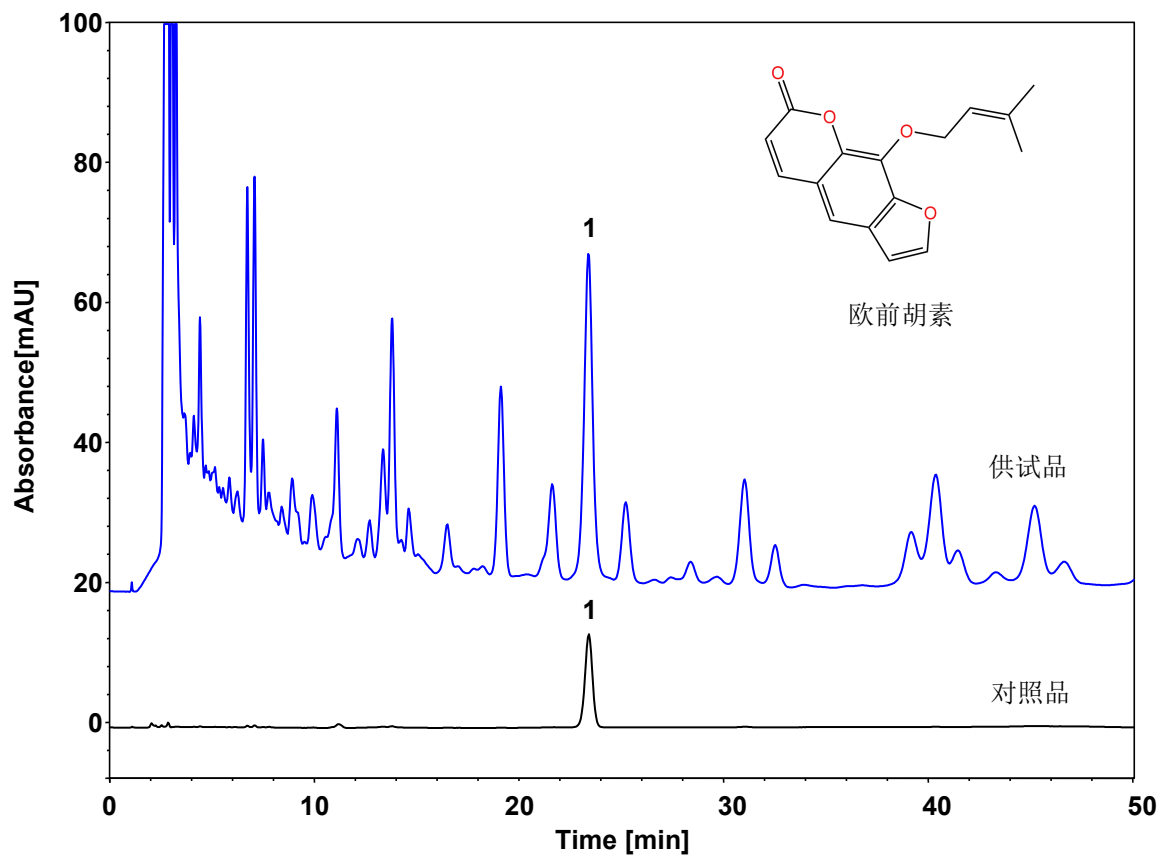

Column: ChromCore AR C18, 5  $\mu$ m  
Dimension: 4.6  $\times$  250 mm  
Mobile Phase: 58/42 v/v 甲醇/水  
Flow Rate: 1.0 mL/min  
Temperature: 35  $^{\circ}$ C  
Injection: 10  $\mu$ L  
Detection: UV 300 nm  
Instrument: Waters 2695  
Sample: 花粉祛痒止痛酊  
Peak: 1. 欧前胡素

组分	保留时间 (min)	理论塔板数 (USP)	拖尾因子 (USP)	分离度 (USP)
欧前胡素	23.395	15460	0.98	2.37