


大黄 (中国药典)

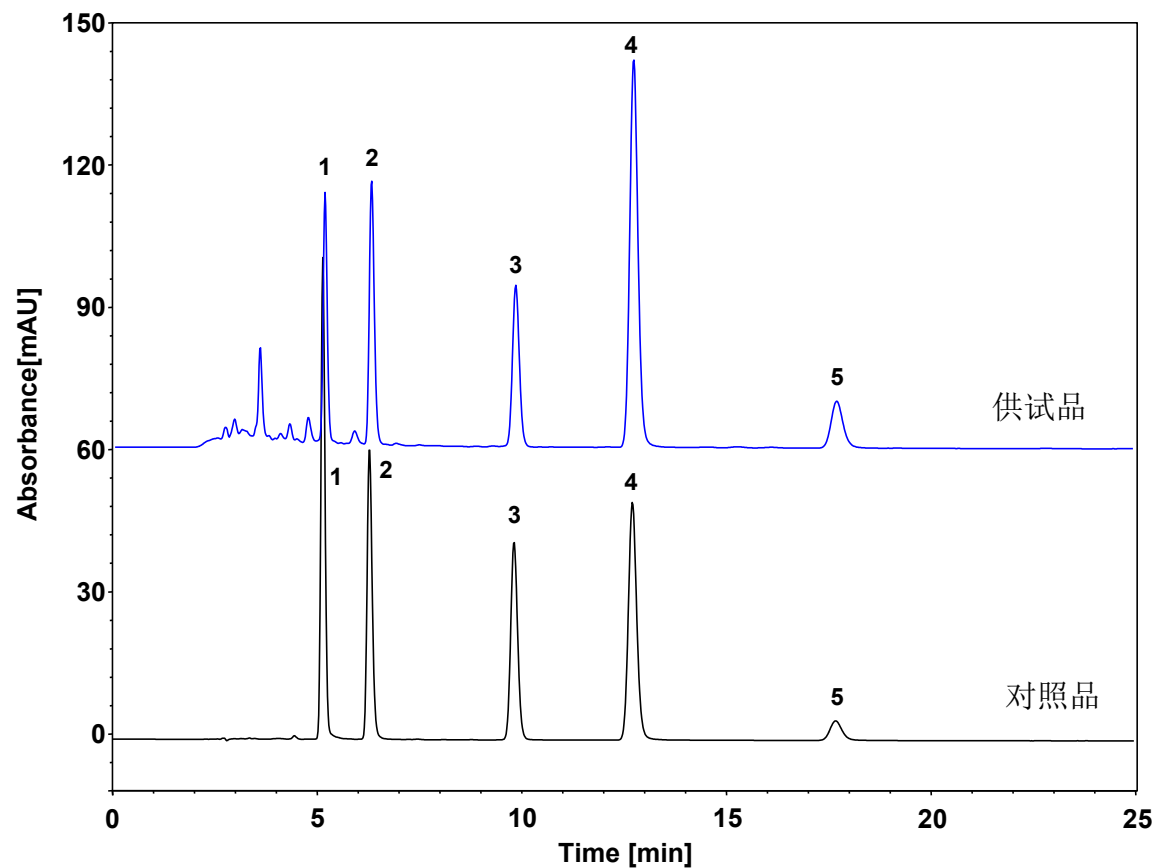

Column: ChromCore 120 C18, 5 μ m
Dimension: 4.6 \times 250 mm
Mobile Phase: 85/15 v/v 甲醇/0.1%磷酸水
Flow Rate: 1.0 mL/min
Temperature: 25 $^{\circ}$ C
Injection: 10 μ L
Detection: UV 254 nm
Sample: 大黄
Peaks:
1. 芦荟大黄素
2. 大黄酸
3. 大黄素
4. 大黄酚
5. 大黄素甲醚