


# 抗体-效价分析

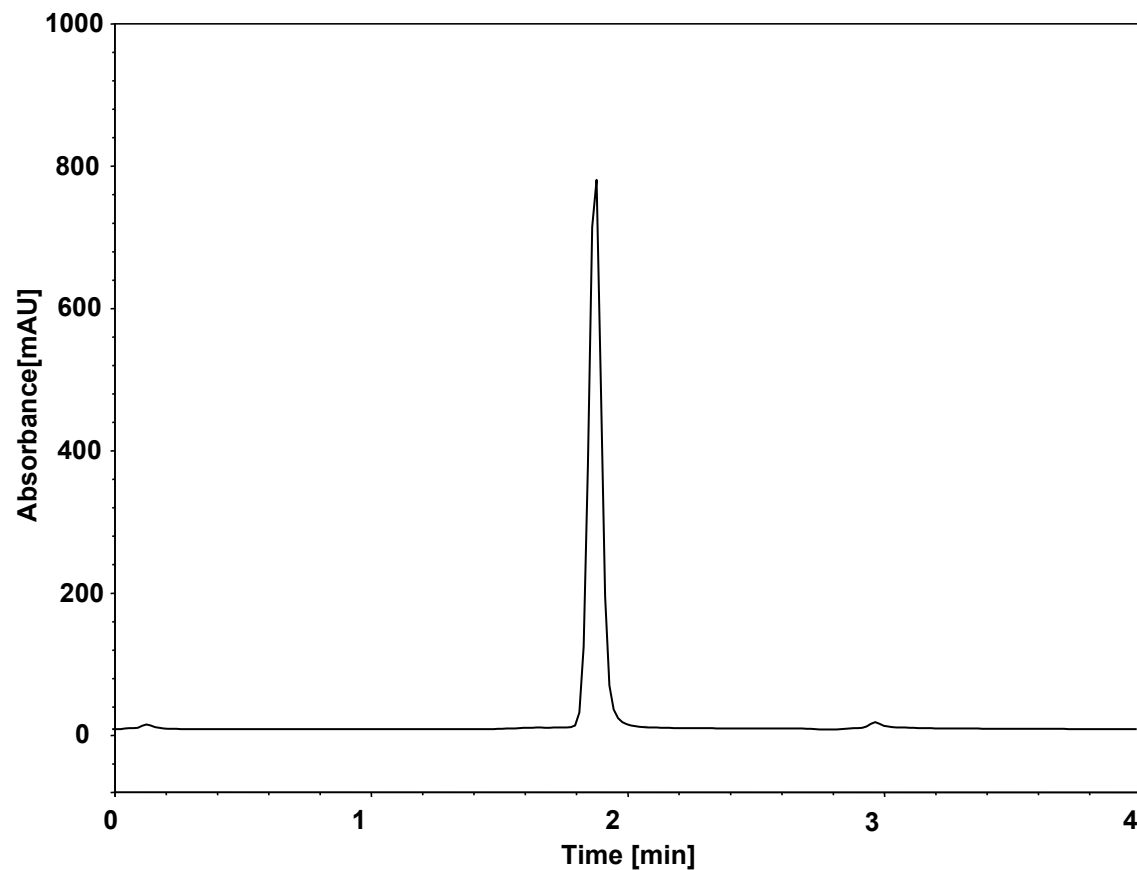

Column: BioCore Protein A, 15  $\mu$ m  
Dimension: 2.1  $\times$  50 mm  
Mobile Phase: A) 150 mM NaCl in 50 mM phosphate buffer, pH7.0  
B) 150 mM NaCl in 50 mM phosphate buffer, pH2.5  
Gradient: 

t (min)	%A	%B
0	100	0
1	100	0
1.01	0	100
2.2	0	100
2.21	100	0
4	100	0

  
Flow Rate: 2.0 mL/min  
Temperature: 25  $^{\circ}$ C  
Injection: 20  $\mu$ L  
Detection: UV 280 nm  
Sample: 抗体